




<b>Technologies</b>    	Product Classification: <b>PLAY &gt; Progressive</b>	<b>Marking:</b> 
	<b>Minimum Fitting Heights:</b> 19mm	

### The Lens Design

Every golfer has his day... you know the day we're talking about. The sun was shining, the birds were chirping, your swing was powerful and exact, and the ball did everything you willed it to... you played the best game of your life. Now with Shamir **Golf™**, every day out on the green can be your best day!

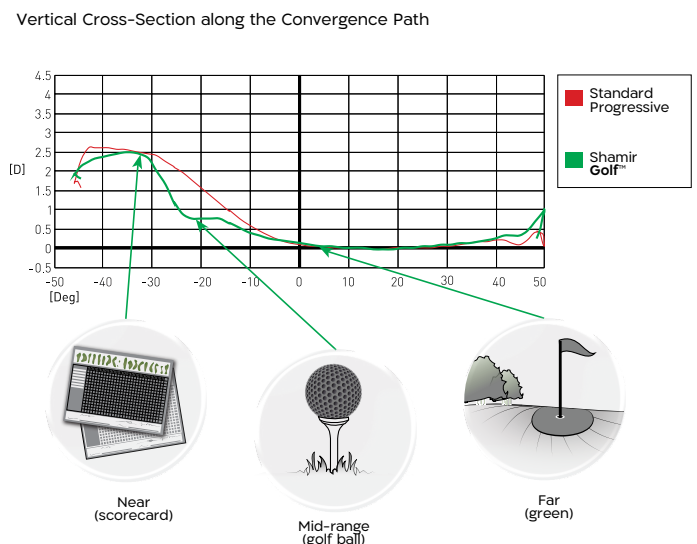
### Designed For Your Visual Sweet Spots.

Shamir **Golf™** is a **Freeform®** lens designed to ensure all your patients experience their perfect game. Unlike standard progressives, the three vision zones are specifically tailored to the needs of golfers to maximize their viewing: the scorecard in their hand, the ball at their feet and the green in the distance. Shamir **Golf™** provides sharp focus without distortion in these three crucial areas. And before we forget, it's the first Golf lens designed with a stable mid-range zone to ensure the ball on the ground can be seen clearly so your patients will enjoy more accuracy whether driving from the tee or using irons on the fairway.

### Who It's For

Avid golfers who need clear vision that doesn't disrupt their focus, or their game.

### One? Two? THREE Enhanced Vision Zones!



IntelliCorridor™ technology makes these enhanced vision zones possible! Learn more about IntelliCorridor™ below.

**As-Worn Technology™***Perfect Vision For Any Frame Choice*

With As-Worn Technology™, we use three different measurements to create a lens that's even better suited for the patient; with vertex distance, pantoscopic tilt and panoramic angle we can fine-tune the lens unlike ever before.

**EyePoint Technology™***The Design Inside*

EyePoint Technology® simulates the human eye in every angle, prescription, and field of vision. These techniques enable Shamir to create the most sophisticated progressive lens surfaces based upon thousands of points of data. It's this "design inside" that makes Shamir progressive lenses the most advanced in the world.

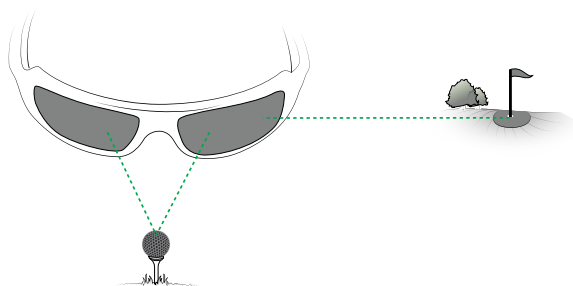
**IntelliCorridor™***Unique Power Profile For Clearer Vision*

Makes it possible to customize lenses specifically to the patient's vision requirements. Optimizes the power in specific areas used for everyday and outdoor tasks, reducing the need to move head up and down.

**Why It's Important**

A golfer's peripheral vision is extremely important when they are lining up to take a swing. They need to be able to focus on the ball at their feet and at the same time see where they want the ball to go in the distance. This is where they utilize their peripheral zone and need it to be crystal clear. Lucky for them, Shamir Golf™ is designed specifically for wraparound frames and has a clear peripheral vision zone so your patients can keep on improving their game.

Thanks to EyePoint Technology™, and As-Worn™ measurements, crystal clear vision on the course is made possible.

**The Secret is in the Peripheral Vision.**

Take it from the experts who designed the first distortion-free solution for wrap frames, when it comes to providing those patients who enjoy spending time on the green with a progressive solution, Shamir Golf™ will give them the vision they need.

**Measurement Information****As-Worn Technology™**

To get the maximum benefit from As-Worn Technology™, the actual measurements of the wearer should be included. (Vertex Distance - Fitted, Pantoscopic Tilt, Panoramic Angle). **The patient will benefit greatly from using their personal measurements, experiencing more precise and comfortable vision.**

If the As-Worn™ measurements are not included with the Shamir Golf™ order, the Shamir Prescripator Software® will use defaults. These defaults are based on averages.

**Default Measurements:**

Vertex Distance: 13 mm  
Pantoscopic Tilt: 9°  
Panoramic Angle: 15°

**"As-Worn" measurements:**

Vertex Distance: 5 to 30 mm  
Pantoscopic Tilt: -10° to 20°  
Panoramic Angle: 13° to 30°



/ShamirInsight



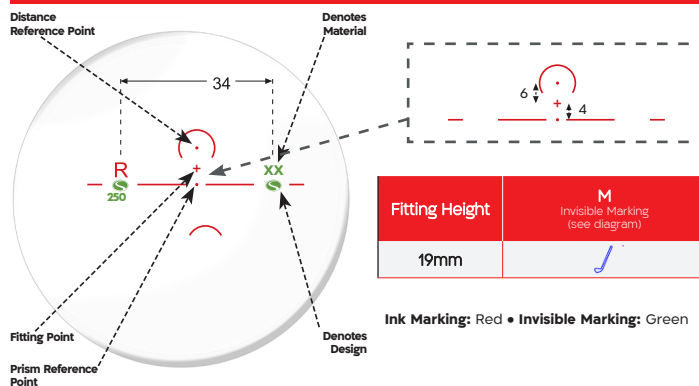
@Shamir\_Insight



thevitaminsee.com



SHAMIRLENS.COM 877.514.8330

**Technical Information****Availability:**

Material	Power Range	Add Power	Cylinder
DLC™ Trivex™ (C,T,V,T7,P,TX,BZ)	-5.00 to +4.00	1.00 to 3.50	to -6.50
1.56 (C,P)	-5.50 to +5.25	1.00 to 3.50	to -6.75
1.6 SuperLite™ (C,T7,P,TX)	-6.00 to +5.25	1.00 to 3.50	to -7.00
Polycarbonate (C,T,V,T7,P,TX,D,BZ)	-6.00 to +5.25	1.00 to 3.50	to -7.00
NXT (PH)	-5.00 to +4.00	1.00 to 3.50	to -6.50

**LEGEND:**

- C - Clear
- TV - Transitions® Vantage™
- T7 - Transitions® Signature7™
- P - Polarized
- TX - Transitions® XTRActive™
- B - BLUtech (Indoor & Outdoor)
- D - Drivewear®
- N - NXT™
- BZ - Blue Zero®
- S - SunSync®
- PH - Photochromic

POWER RANGE IS DETERMINED BY SHAMIR CERTIFIED BLANK RANGE PRISM TO 10 DIOPTERS IN EACH QUADRANT

**⚠ WARNING:** Polycarbonate lenses can expose you to chemicals including Bisphenol A (BPA) and 1,74 index lenses can expose you to Methimazole, which are known to the State of California to cause birth defects or other reproductive harm. For more information go to: [www.P65Warnings.ca.gov](http://www.P65Warnings.ca.gov).

**Features & Benefits**

- Specialty lens designed specifically for the golf enthusiast
- Designed to ensure sharp vision at 3 crucial focal points
- Utilizes As-Worn Technology™ and EyePoint Technology®
- Provides the player with a more comfortable golfing experience
- Clear peripheral viewing
- Personalized to the patient
- Designed specifically for wrap frames

**How To Fit & Order Shamir Golf™****Provide the following information:**

- Patient Rx
- Fitting Height
- Monocular Distance PD
- Panoramic Angle
- Pantoscopic Tilt
- Fitted Vertex Distance
- Frame Measurements A/B/DBL&ED

If the As-Worn™ measurements are not included with the order the Shamir Prescripator Software® will use defaults. Please see the measurements section for more details.

\*It is strongly recommended that you provide all of the above mentioned measurements.